


E. coli



Purple, with purple halo

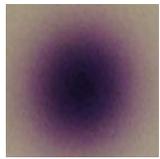

Purple, no halo



Purple with pink halo



Blue with purple or pink halo



Blue or dark blue, no halo



Dark blue with teal halo



Dark blue with blue halo

Not E. coli



White




Pink, no halo




Pink with pink halo



Pinpoints*
(If after incubation period)

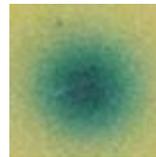

Teal green, no halo

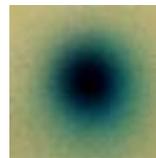

Teal with teal halo



Red

*Do not count pinpoints if the plate is dominated by larger colonies. Pinpoints may be counted if they make up >50% of colonies. If possible, incubate a few additional hours to see if colonies will grow larger